ToothWave

H7001

- EN USER MANUAL
- FR MODE D'EMPLOI
- ES MANUAL DEL USUARIO

Patent protected
technology by

TOOTH
WAVE™

SCIENCE DRIVEN ORAL CARE

DentalRF™
Technology

Silk'n ToothWave

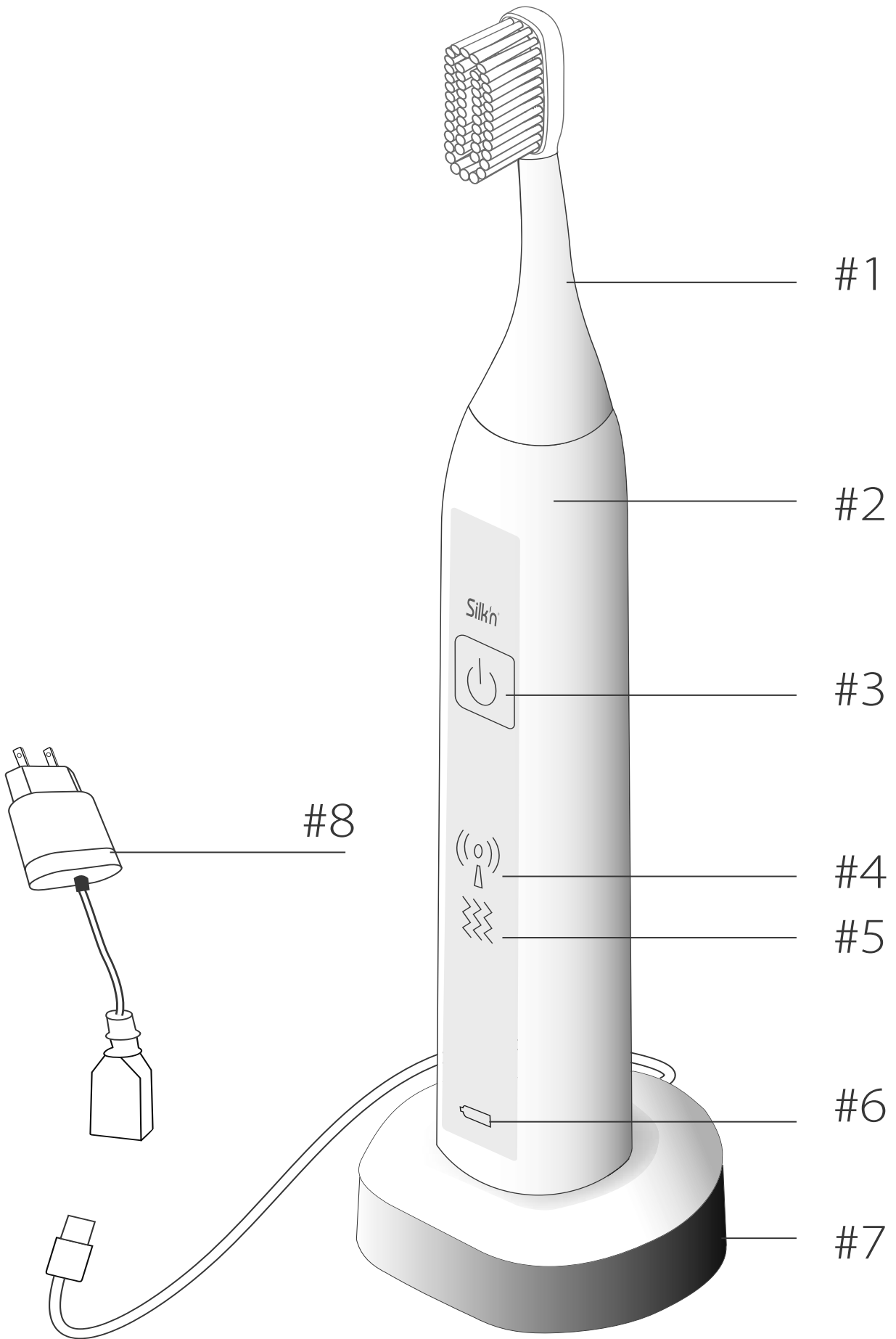
Table of contents

1.	Warnings and precautions	2
2.	Contraindications	5
3.	Get to know the device.....	6
3.1	Indications for use	6
3.2	Clinical summary data	6
3.3	Package contents.....	6
3.4	Device controls & settings.....	7
3.5	The brush head	8
3.6	Device charging.....	9
4.	Treatment plan.....	9
5.	What to expect with the device	10
6.	Get started!.....	10
7.	Maintenance.....	11
7.1	Routine maintenance.....	11
7.2	Maintenance of the RF contact.....	11
7.3	Cleaning the charging cradle	12
7.4	Replacing the brush head	12
8.	Parts & consumables	13
9.	Troubleshooting	13
10.	Customer service	14
11.	Warranty	14
12.	Disposal	14
12.1	Battery removal.....	14
12.2	Packaging.....	17
12.3	Device	17
13.	Specifications	17
14.	Signs & symbols	19
15.	Labeling	20
16.	Output power curves	20
17.	EMC Technical results	21
18.	FCC Declaration of Conformity	22
19.	Legal notice	23

1. Warnings and precautions

- Read the manual in full before use, and refer to it if you are in any doubt about how to use the device safely.
- The device is for home-use. The handle may be used by multiple persons, but for reasons of hygiene, each user must use their own brush head (see section [7.4 Replacing the brush head](#)).
- Use of this device is not a substitute for regular visits to the dentist for routine clinical care.
- The device is not intended for use on multiple patients in a dental practice or institution. Do not use for any purpose other than as described in this manual.
- Always ensure that your hands are dry before connecting the charging cradle to the mains supply.
- Do not place or store the USB wall adapter where it can fall, or be pulled into any liquid.
- Do not use the USB wall adapter if it has fallen into liquid.
- For the purpose of recharging the battery, only use the USB wall adapter provided with this device.

- If the charging cradle malfunctions, unplug immediately.
- Never use the charging cradle or the USB wall adapter if they are damaged or aren't working properly.
- Keep the device, all parts and accessories away from heated surfaces.
- Do not attempt to modify, open or repair the device. You may expose yourself to an electric shock.
- Do not use if the parts or accessories are not supplied by the manufacturer as defined in this manual, or if they are damaged in any way, do not appear to be working properly, or if you see or smell smoke. In such cases, stop using the device, and contact customer service.
- The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.
- Do not clean any of the parts in the dishwasher.
- Keep the device out of the reach of children or pets.
- This device contains batteries that are non-replaceable.

- ⦿ Do not dispose of the device or the battery with household waste! The law requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices.
- ⦿ Dispose of the battery safely.
- ⦿ If the product reaches its end of life, or you no longer wish to use it, the battery must be removed before disposal (see section 12. Disposal).

Precautions

⚠ Consult your physician before use if you:

- ⦿ have concurrent conditions, such as heart disorders, seizures, uncontrolled high blood pressure, liver or kidney diseases.
- ⦿ have poorly controlled hormone disorders, such as diabetes or thyroid abnormalities.

⚠ Consult your dentist before use if you:

- ⦿ have had oral or gum surgery in the last two months.
- ⦿ have severely injured gums.

 **Additional precautions:**

- Do not bite the brush head.
- Do not touch or apply pressure to the brush head when you use, clean or maintain the device.
- The safety and effectiveness of the device for uses other than those indicated are not known.

2. Contraindications

This device including all its parts must not be used by the following people: persons below the age of 18, those with limited physical, sensory or psychological capacities, those lacking experience or knowledge in how to use the device in a safe way, or those who do not understand the hazards involved.

 **Do not use the device if you:**

- have a pacemaker, internal defibrillator, another active implanted device, or if you have any medical concerns about using this device.
- have a history of oral cancer, throat or oropharyngeal cancer, or have any other type of cancer, or have pre-malignant lesions.
- have an impaired immune system due to illnesses, or you take medication that impairs the body's ability to fight infection.
- are pregnant or nursing.

3. Get to know the device

EN

3.1 Indications for use

ToothWave™ is a powered radiofrequency toothbrush intended to promote good oral hygiene, including reduction of plaque and the prevention and treatment of gingivitis. ToothWave™ is intended for over-the-counter use.

3.2 Clinical summary data

In a 6 week clinical study, the control group of 42 subjects using a powered ADA Accepted toothbrush (31,000 strokes/minute) showed a 15% reduction in plaque, a 31% reduction in gum inflammation and a 2.8% reduction in gum bleeding, while the 45 subjects using the ToothWave™ (33,000 strokes/minute) showed a 25% reduction in plaque, a 45% reduction in gum inflammation and a 41% reduction in gum bleeding.

In the same study where patients began with a moderate amount of visible tartar on the tongue side of their lower anterior teeth, additional analysis showed a 14% accumulation of tartar in the control group and a 2% accumulation of tartar in the ToothWave™ users.

3.3 Package contents

Refer to the diagram on page (i) at the front of the manual.

- 2x Brush heads (#1)
- 1x Handle (#2)
- 1x Charging cradle with USB cable + transparent base (#7)
- 1x USB wall adapter (#8)
- 1x Warranty booklet
- 1x User manual)

3.4 Device controls & settings


Refer to the diagram on page (i) at the front of the manual.

Control	Function / Mode
Control button (#3) 	Long press: turns the device on or off Short press: each press moves to the next mode
DentalRF mode indicator light (#4) 	DentalRF is active for all modes 
Vibration mode indicator light (#5) 	Mode 1: low vibration 
	Mode 2: medium vibration 
	Mode 3: high vibration Recommended mode 
	Mode 4: No vibration 
Battery indicator light (#6) 	Green blinks: battery is charging No light: battery is charged Orange blinks: battery is low Orange constant: battery is empty

i When the battery is fully charged, the light will be constant green for 30 seconds, then it will go out.

3.5 The brush head

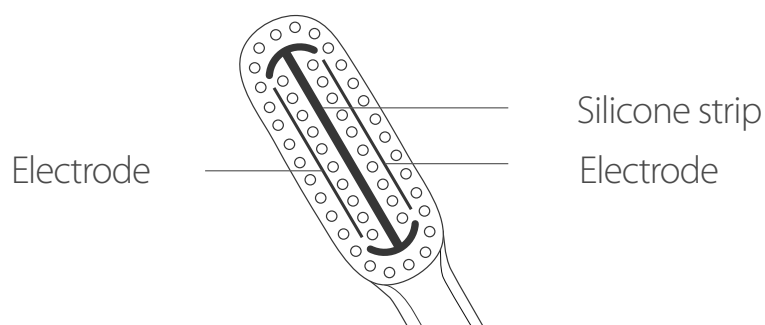The brush head **(#1)** is an innovative and special feature of the device. It contains bristles, 2 DentalRF electrodes and a silicone strip that lies between the electrodes. The silicone strip enables the DentalRF current to flow across the teeth. It also ensures that the toothpaste and the DentalRF energy interact effectively. For this reason, we recommend that the bristles are always in good contact with the teeth.

- ⚠ Always check that the brush head **(#1)** is not damaged. Inspect the brush head **(#1)** for damage, excessive wear, and/or biofilm accumulation prior to each use, and replace the brush head **(#1)** or handle **(#2)** if there are signs of damage.
- ⚠ If you are in any doubt about its condition, replace with a new one. Naturally, contact Customer Service at any time for specific queries or questions.

Signs of damage to the brush head **(#1)** include but are not limited to:

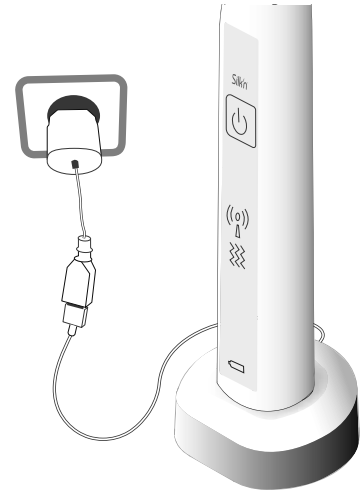
- a significant change in shape or appearance;
- excessive wear showing on the bristles;
- a bend or tears in the silicone barrier.

Do not use the brush head **(#1)** if the bristles, silicone strip, or electrodes become loose, bent, damaged, or crushed. Damaged electrodes, bristles or a damaged silicone strip may break off while brushing. If it is damaged, replace with a new brush head **(#1)**. Refer to section [7.4 Replacing the brush head](#).


3.6 Device charging

- Check that the battery is not low before each use.

4. Treatment plan

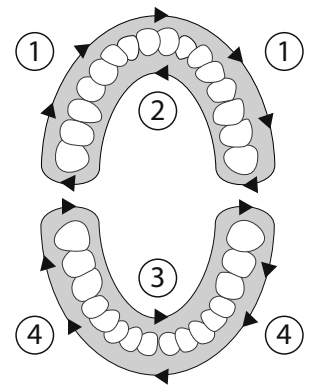- Use for 2 minutes, morning and evening.
 - Clean between the teeth and gently massage gums.
 - Brush each section (1,2,3,4) for 30 seconds.
- i** The device is suitable for use with dental implants.

For details, refer to section [6. Get started!](#).

5. What to expect with the device

During or after using the device you may:

- ⦿ Feel pain, if you have cavities or any teeth condition that requires dental treatment. If you do feel pain, consult your dentist.
- i** Stop using the device if you experience excessive bleeding after use, or if bleeding continues after one week of regular use.

6. Get started!

1. Rinse the bristles under lukewarm water.
2. For best results, apply a generous amount of toothpaste. You can use any commercial toothpaste.
3. Hold the handle **(#2)** firmly, and turn on the device. Short press the control button **(#3)** as many times as needed to select your treatment mode.
4. Gently press the bristles against the teeth, and move the brush in circular movements, slowly and systematically from tooth to tooth, gently massaging the gums. Make sure the bristles are always in good contact with your teeth.

i The device comes with soft bristles to protect your teeth and gums. Brush gently, and apply only light pressure when brushing. Do not scrub.

5. Brush all surfaces of the teeth for 2 minutes. You will feel a slight vibration every 30 seconds (1 vibration at 30 seconds, 2 vibrations at 60 seconds and so on). After the vibration(s), you can move to another section of the mouth (refer to section [4. Treatment plan](#)). After 2 minutes, the device will shut down automatically.

- Clean the device with warm water, and store in open air whilst it is wet.

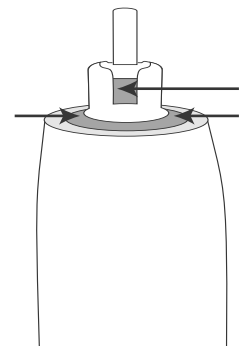
7. Maintenance

7.1 Routine maintenance

- After use, wash the brush head under running water to remove any toothpaste or debris.
- Do not wash with soap or chemicals.
- If debris or toothpaste accumulates on the handle (#2), wash under running water and dry with a clean dry cloth.
- Do not try to remove any stuck debris from the brush head (barrier and/or electrodes). In such cases, replace the brush head (#1). Refer to section [7.4 Replacing the brush head](#).


7.2 Maintenance of the RF contact

- Once a month, remove the brush head (#1) from the handle (#2). Clean the groove in the handle (#2) to remove any toothpaste or debris that may have accumulated. Use a cotton bud, rinse and dry with a clean dry cloth.
- Remove debris and excess toothpaste from the brush head (#1) and handle (#2). Clean the underside of the brush head (#1) using a cotton bud.
- If another person is going to use their brush head (#1), be sure to rinse and wipe down the handle (#2), including the connection between the brush head (#1) and the handle (#2) before another person uses the device.

7.3 Cleaning the charging cradle

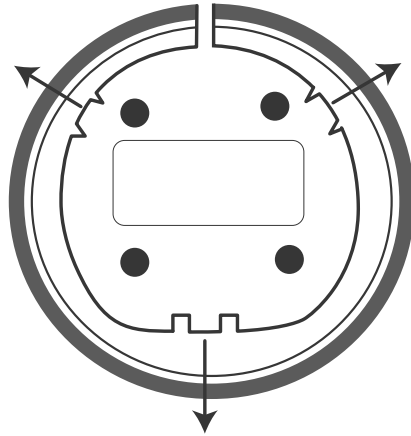1. First disconnect the charging cradle (**#7**) from the mains supply and remove the toothbrush from the charging cradle (**#7**).
2. Next, remove the **transparent base**. To do this, turn the charging cradle (**#7**) upside down. Remove the transparent base by pressing the clips outwards.


3. Clean both parts with a damp cloth and dry them. Place the charging cradle (**#7**) back, and press down firmly. You will hear it click into place.

7.4 Replacing the brush head

Replace the brush head (**#1**) when any of these apply:

- another person wishes to use the device, or
- you have used the same brush head (**#1**) for 3 months, or
- the bristles are damaged or worn.

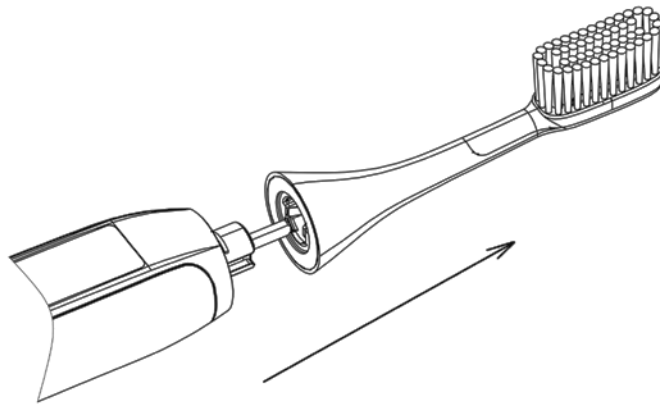To replace, simply pull out the brush head **(#1)**, remove the old one and replace with a new one.

8. Parts & consumables

You can purchase replacements from where you purchased the device, or online. Refer to the Customer Service section for details.

9. Troubleshooting

Do not attempt to open or repair your device. You may expose yourself to dangerous electrical components and cause serious harm to yourself. It will also void your warranty.

Problem	Checks
My device does not start.	Make sure the battery is not low. Refer to section 3.6 Device charging .
The device is charged but it does not work.	Press the control button (#3) to restart the device. If the problem persists, contact customer service.
All indicator lights (#4, #5, #6) blink.	There is a system error. Press the control button (#3) to restart the device. If the problem persists, contact customer service.

i For problems or queries, please contact Silk'n Customer Service.

10. Customer service

EN See our website for more information: www.silkn.com. This manual is also available as a PDF download from: www.silkn.com/manuals. Contact Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

✉ contact@silkn.com

☎ 1-877-367-4556 (toll-free in the US)

11. Warranty

Please refer to the separate product warranty booklet for complete information.

i The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.

12. Disposal

Be sure to follow your country's laws for the disposal of electrical items. Contact your local council for information on recycling points and/or disposal services for lithium batteries and electrical devices.


12.1 Battery removal

At the product's end of life, you are required to remove the battery before disposal. To do this, the device must be opened up. The device will no longer function after this procedure, nor can it be repaired. Any warranty will be void when the device casing is opened and the battery is removed. Home Skinovations Ltd is not liable for any physical or material damages as a result of performing this procedure.


Refer to the diagrams at the end of this section for steps 5-8.

1. Remove the device from the charging cradle **(#7)**.
2. Make sure that the battery is completely empty. Turn on the device to fully drain the battery.
3. Remove the brush head **(#1)** from the handle **(#2)**.
4. As a protective measure, wear safety goggles for the rest of the process.
5. Open the casing: Put the flat-head screwdriver into the slot at the bottom of the handle **(#2)** and gently swing the bottom cover gently until you see a gap - see diagram step 5.
6. Place the flat-head screwdriver into the gap and pry the bottom section out of the handle **(#2)** - see diagram step 6.
7. Remove the inner section: Place the handle **(#2)** upside down on a sturdy table covered with a towel. Use the flat-head screwdriver to move the clips inwards. As the clips release, push down to release the inner part of the handle **(#2)** - see diagram step 7.
8. Remove the battery:
 - ⚠ If the battery is damaged or leaking, wear rubber gloves before touching it.
 - ⚠ Make sure the exposed red and black wires do not touch each other. Cut each wire separately.
 - ⚠ If the red and black wires touch each other, separate them immediately.
 - With cutters, cut the red wire of the battery close to the printed circuit board (PCB). Seal the exposed end (battery side) with insulating tape - see diagram step 8.

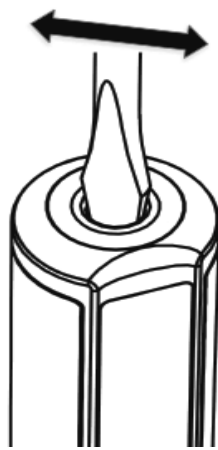
- Again with cutters, cut the black wire of the battery close to the PCB. Seal the exposed end (battery side) with insulating tape - see diagram step 8.
 - Carefully remove the battery from the casing.
9. Place the battery in a plastic bag: This prevents the battery from coming into contact with other items.
 10. Dispose of the battery and the device as advised by your local council.

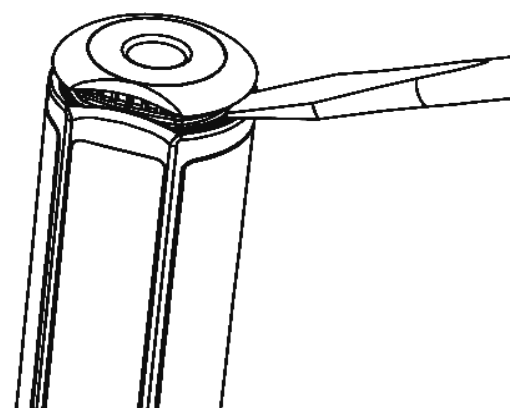
Step 5

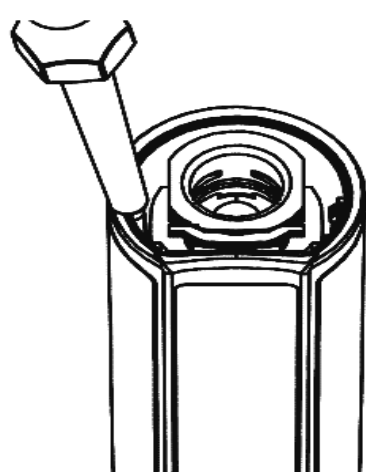Step 6

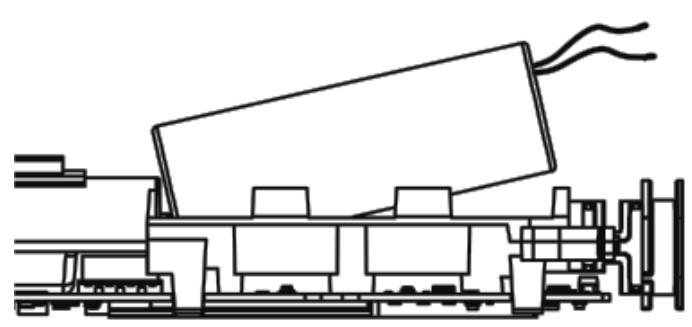Step 7

Step 8

12.2 Packaging

Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

12.3 Device

Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. This ensures that the device is recycled professionally, and prevents the release of harmful substances* into the environment. Remove the battery and dispose of it safely before disposing of the device.

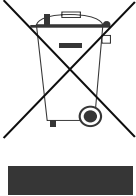
* Products labelled with Cd = cadmium, Hg = mercury, Pb = lead.

13. Specifications

Model no.	H7001
Technology	DentalRF™ and vibration
Vibration	0; 275Hz; 300Hz or 400Hz (± 5%)
Radio frequency	3MHz ±0.3MHz; up to 3W
Package size	(W)165 (H)227 (D)80 [mm]
System weight	115g
Transport & storage between uses and storage condition	Temperature: - 40 to 70°C Relative humidity : 10 to 90%rH Atmospheric pressure: 500 to 1060hPa
Operating conditions	Temperature: 5 to 40°C Relative humidity: 15 to 90%rH Atmospheric pressure: 700 to 1060hPa

Time from minimum storage temperature between uses to minimum operating temperature with ambient temperature of 20°C	30 minutes
Time from maximum storage temperature between uses to maximum operating temperature with ambient temperature of 20°C	30 minutes
Protection against electric shock	Class I medical equipment type BF applied part
Adapter	YH-S06U0500600 (USA)
Input: Output:	100 – 240V; 50/60Hz; 0.2A 5.0Vdc; 0.6A
Mode of operation	Non-continuous
Waterproof rating:	
USB adapter (IPX4)	Protection assured against water splashing
Charging cradle (IPX7)	Protection assured against water immersion at a depth of up to 1 meter, for a maximum of 30 minutes.
Handle (IP67)	Totally protects against dust and protects in water up to 3 feet (1 meter).
The device was tested and complies with	IEC 60601-1 (Safety) IEC 60601-2-2 (RF Safety) IEC 60601-1-2 (EMC) IEC 60601-1-11 (Environment)
Service life of the equipment	5 years

14. Signs & symbols

	<p>Important information about using or maintaining the device.</p>
<p>IP67</p>	<p>Protected from immersion in water with a depth of up to 1 meter (or 3.3 feet) for up to 30 minutes.</p>
	<p>Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.</p>
<p>Li-ion</p>	<p>Lithium-ion battery inside.</p>
	<p>Degree of protection against electric shock: type BF applied part.</p>
	<p>Follow the operating instructions.</p>
	<p>The built-in rechargeable battery contains substances that may pollute the environment. Always remove the battery before you discard at an official collection point for batteries.</p>
	<p>This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.</p>
	<p>Manufacturer</p>
	<p>Recycle the packaging.</p>
	<p>Worldwide warranty.</p>

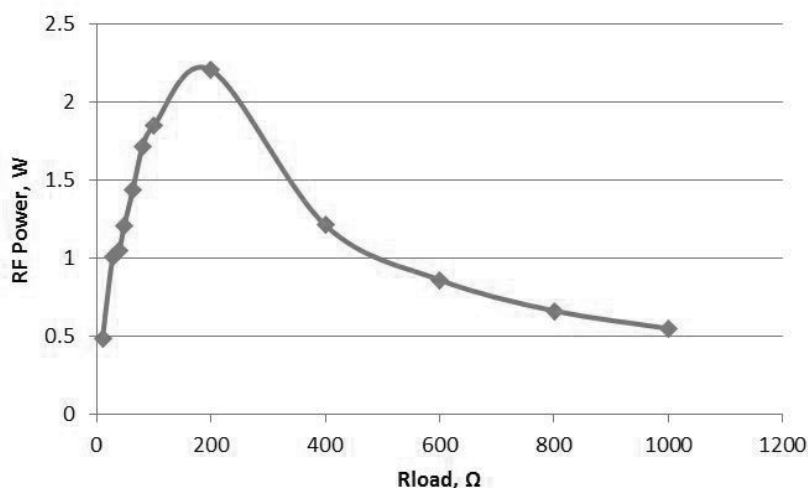
15. Labeling

This device is not suitable for use in the presence of flammable anesthetic mixtures with air, oxygen or nitrous oxide.

16. Output power curves

The load (impedance) is the resistance of the toothpaste and saliva medium through which the RF current flows. The RF power is the rate per second at which the RF energy is transferred by the RF generator into the impedance. During brushing, the toothpaste density is changed, as the toothpaste is foamed and mixed with saliva, this change in density changes the impedance respectively. The curve shows the device RF power (for all modes of operation) versus the toothpaste and saliva medium impedance.

17. EMC Technical results

EMC requirements

IEC 60601-1-2 Edition 4.0 (2014)

Environment of intended uses:

Home Healthcare Facility Environment

Summary of Test Results

Test	Standard	Class/Severity level	Test result
Emission			
Conducted emission Freq. range: 150 kHz – 30 MHz	FCC Part 15 Subpart B	Class B on 120 Vac mains	Complies
Radiated emission Freq. range: 30 MHz – 1 GHz		Class B	Complies
Voltage changes, Voltage fluctuations and Flicker test	EN 61000-3-3	AC mains	Complies
Immunity (IEC 60601-1-2 section 8)			
Immunity from Electrostatic discharge (ESD)	IEC 61000-4-2	8 kV contact discharges & 15 kV air discharges	Complies
Immunity from radiated electromagnetic fields	IEC 61000-4-3	3.0; 10.0 V/m; 80 MHz 2.7 GHz, 80% AM, 1 kHz	Complies
Immunity from Proximity field from wireless communications equipment	IEC 61000-4-3	List of frequencies, from 9 V/m up to 28 V/m, PM (18 Hz or 217 Hz), FM 1 kHz	Complies
Immunity from power frequency magnetic field	IEC 61000-4-8	30 A/m @ 50 Hz	Complies

18. FCC Declaration of Conformity

EN

Home Skinovations Ltd.
Tavor building,
Shaar Yokneam, POB 533
Yokneam 2069206, Israel
www.silkn.com

November 26, 2020

Declaration of Conformity 47 CFR § 2.1077 Compliance Information

Model: ToothWave H7001

Responsible Party U.S. Contact Information:

Yair Malca
Inmode U.S. office
17 Hughes Irvine, CA 92618
USA

Email: yair.malca@inmodemd.com

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by Home Skinovations Ltd. could void the user's authority to operate the ToothWave device.

Signature:

Date: November 26, 2020

Name: Ben Zion Levi

Title: co-CEO

19. Legal notice

Copyright 2020 © Home Skinovations Ltd. All rights reserved.

Home Skinovations Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification. Home Skinovations Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice. Silk'n and the Silk'n logo are registered trademarks of Home Skinovations Ltd, Tavor Building, Shaar Yokneam P.O. Box 533, Yokneam 2069206, ISRAEL.

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Copyright 2020 © Home Skinovations Ltd. Tous droits réservés.
Copyright 2020 © Home Skinovations Ltd. Todos los derechos reservados.